

CHARACTERIZATION OF THE INTEGRATED VIRUS DETECTION SYSTEM (IVDS): A NOVEL INSTRUMENT FOR VIRUS IDENTIFICATION AND ITS USE IN THE INITIAL CHARACTERIZATION OF A NEW MARINE VIRUS

D.A. Kuzmanovic*

Geo-Centers, Inc., Aberdeen Proving Ground, MD 21040

I. Elashvili and C. H. Wick

Edgewood Chemical Biological Center, APG, MD21010-5424

S. Krueger

NIST Center for Neutron Research, Gaithersburg, MD 20899-856

ABSTRACT

Detection and quantification of unknown or novel virus samples remains one of the great technical challenges in biological sciences. Here we report the characterization of an instrument, The Integrated Virus Detection System (IVDS) as an accurate and rapid means to size and quantify viral unknowns.

1. INTRODUCTION

Typically, virus size is determined by electron microscopy and virus concentration is measured by infectivity or optical density (O.D.) if the virus host and other physical features of the virus are known. Recent outbreaks of previously undescribed viral infections and the announcement of the creation of recombinant viruses have highlighted the need for a single instrument that can detect and identify viruses in solution about which nothing is known. Here, the Integrated Virus Detection System (IVDS), the only virus identification system designed for virus quantification and sizing, is analyzed in terms of the accuracy and reliability of its measurements (Wick, 1999). Specifically, to test the quantification ability of the IVDS instrument for virus sizing and counting two approaches were employed: 1) The accuracy of sizing by IVDS was determined using synthetic and biological sizing standards. 2) Evaluation of the concentration determination ability of IVDS was determined using the MS2 bacteriophage as a biomarker and Small Angle Neutron Scattering (SANS).

2. RESULTS AND DISCUSSION

2.1 SIZING ANALYSIS BY IVDS

The accuracy of sizing by IVDS was determined using synthetic and biological standards. Specifically, National Institute of Standards and

Technology (NIST)-traceable size standard latex microspheres (data not shown) and the MS2 bacteriophage were both measured by IVDS. See Figures 1.

Figure 1.

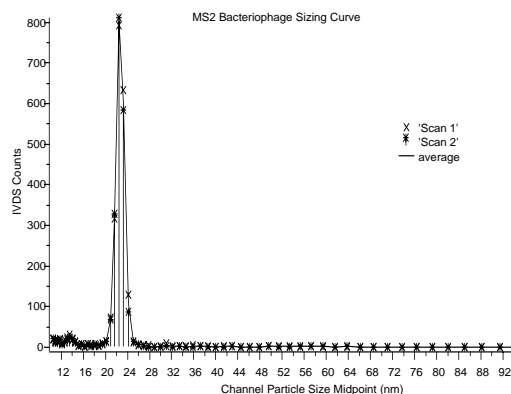


Table 1: Summary of MS2 Bacteriophage Particle Size Determinations

Sizing Method	Size (nm)	References
IVDS	23	This study
TEM	25	Overby, 1966
TEM	24	Sugiyama 1967
TEM	26	Straus 1963

Taken together, this data shows IVDS sizing measurements for the MS2 virus and NIST-traceable size standard latex microspheres are 23 and 43nm, respectively. This is in good agreement with the reported values above.

2.2 IVDS STANDARD CURVE

Also since the NIST-traceable microspheres are extremely well characterized, with respect to size as well as concentration, they were also used to generate a concentration standard curve for IVDS. See Figure 2.

Report Documentation Page			Form Approved OMB No. 0704-0188		
Public reporting burden for the collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington VA 22202-4302. Respondents should be aware that notwithstanding any other provision of law, no person shall be subject to a penalty for failing to comply with a collection of information if it does not display a currently valid OMB control number.					
1. REPORT DATE 00 DEC 2004		2. REPORT TYPE N/A		3. DATES COVERED -	
4. TITLE AND SUBTITLE Characterization Of The Integrated Virus Detection System (Ivds): A Novel Instrument For Virus Identification And Its Use In The Initial Characterization Of A New Marine Virus			5a. CONTRACT NUMBER		
			5b. GRANT NUMBER		
			5c. PROGRAM ELEMENT NUMBER		
6. AUTHOR(S)			5d. PROJECT NUMBER		
			5e. TASK NUMBER		
			5f. WORK UNIT NUMBER		
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) Geo-Centers, Inc., Aberdeen Proving Ground, MD 21040; Edgewood Chemical Biological Center, APG, MD21010-5424			8. PERFORMING ORGANIZATION REPORT NUMBER		
9. SPONSORING/MONITORING AGENCY NAME(S) AND ADDRESS(ES)			10. SPONSOR/MONITOR'S ACRONYM(S)		
			11. SPONSOR/MONITOR'S REPORT NUMBER(S)		
12. DISTRIBUTION/AVAILABILITY STATEMENT Approved for public release, distribution unlimited					
13. SUPPLEMENTARY NOTES See also ADM001736, Proceedings for the Army Science Conference (24th) Held on 29 November - 2 December 2005 in Orlando, Florida.					
14. ABSTRACT					
15. SUBJECT TERMS					
16. SECURITY CLASSIFICATION OF:			17. LIMITATION OF ABSTRACT UU	18. NUMBER OF PAGES 2	19a. NAME OF RESPONSIBLE PERSON
a. REPORT unclassified	b. ABSTRACT unclassified	c. THIS PAGE unclassified			

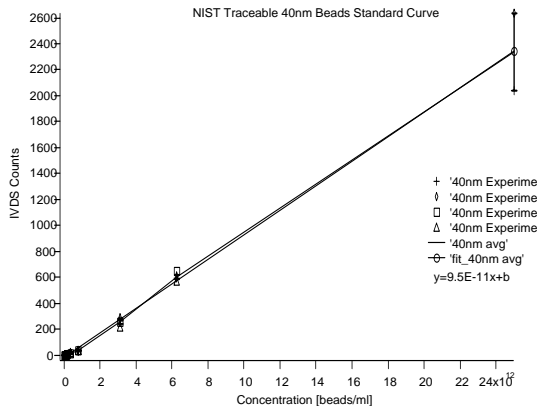


Figure 2.

The NIST microsphere standard concentration curve is linear and sufficient to generate an equation for the conversion of IVDS counts/ml to virus particles/ml.

2.3 EVALUATION OF IVDS CONCENTRATION DETERMINATION ABILITY

Testing of the concentration determination ability of IVDS was determined by measuring the concentration of the bacteriophage MS2 and then using this concentration to determine its molecular weight by SANS.

In general, SANS is a process where a neutron beam is passed through a sample and the resulting scattering pattern reveals information about the average size, shape, and orientation of the sample. Also, when the concentration of the sample is known, then molecular weight can be determined by SANS using absolute scaling (Krueger et. al 1998). Therefore, if IVDS accurately measures virus particle concentration then one should be able to accurately deduce the molecular weight of the virus using concentrations measured by IVDS. Since the molecular weight of MS2 has been determined by a number of different methods, MS2 is an ideal biomarker for use for testing of IVDS.

A stock of MS2 was measure by IVDS to get a concentration in IVDS counts/ml. The data was plotted on the NIST-Traceable concentration curve to convert those numbers into virus particles/ml. The same MS2 stock was also measured by optical density using a conventional spectrophotometry as a control since the use of optical density for virus concentration determination is well established. Finally the MS2 phage stock was subjected to SANS analysis and its molecular weight determined using both IVDS and optical density concentration measurements both methods were

found to be in good agreement yielding resulting in molecular weights of 3.6×10^6 . These results show that IVDS instrument is a rapid and accurate measure of virus sizing and concentration (which unlike other methods) requires no prior information about the unknown virus-containing sample.

3. MATERIALS AND METHODS*

For detailed materials and methods, see (Kuzmanovic et. al. 2003). *Certain commercial materials, instruments, and equipment are identified in this manuscript in order to specify the experimental procedure as completely as possible. In no case does such identification imply a recommendation or endorsement by the NIST nor does it imply that the materials, instruments, or equipment identified is necessarily the best available for the purpose.

REFERENCES

- Krueger, S. (1998). "SANS Provides Unique Information on the Structure and Function of Biological Macromolecules in Solution." Physica B **241**: 1131-1137.
- Kuzmanovic, D. A., I. Elashivili, et al. (2003). "Bacteriophage MS2: Molecular Weight and Spatial Distribution of the Protein and RNA Components by Small-Angle Neutron Scattering and Virus Counting." Structure **11**: 1339-1348.
- Mazzone, H. (1998). CRC Handbook of Viruses: Mass-Molecular Weight Values and Related Properties, CRC Press LLC.
- Strauss, J. H. and R. L. Sinsheimer (1963). "Purification and Properties of Bacteriophage MS2 and of its Ribonucleic Acid." Journal of Molecular Biology **7**: 43-54.
- Sugiyama, T., R. R. Hebert, et al. (1967). "Ribonucleoprotein Complexes formed between Bacteriophage MS2 RNA and MS2 Protein *in vitro*." Journal of Molecular Biology **25**: 455-463.
- Wick, C. H. and P. E. McCubbin (1999). "Characterization of Purified MS2 Bacteriophage by the Physical Counting Methodology used in the Integrated Virus Detection System (IVDS)." Toxicological Methods **9**(4): 245-252.